

[See all 62 Products in Family](#)

## 600 Grooves, 12.5 x 25mm, 400nm Ruled Diffraction Grating



Reflective Ruled Diffraction Gratings

Stock **#41-021** **3 In Stock**

MRP ₹12,107

**Price inclusive of all taxes**

**ADD TO CART**

### Volume Pricing

Qty 1-9	₹12,107 each
Qty 10-24	₹10,896 each
Need More?	<a href="#">Request Quote</a>

### Product Downloads

### General

Reflective Diffraction Grating **Type:**

### Physical & Mechanical Properties

12.5 x 25.0 ±0.5 **Dimensions (mm):**

90 **Clear Aperture (%):**

Ruled Grating	<b>Construction:</b>
Parallel to Short Dimension	<b>Direction of Grooves:</b>
25.00	<b>Length (mm):</b>
9.50 ±0.5	<b>Thickness (mm):</b>
12.50	<b>Width (mm):</b>
±0.5	<b>Alignment of Grooves to Edge (°):</b>

## Optical Properties

600	<b>Groove Density (grooves/mm):</b>
400	<b>Blaze Wavelength (nm):</b>
6.88	<b>Blaze Angle (°):</b>
Bare Aluminum	<b>Coating:</b>
Float Glass	<b>Substrate:</b> <input type="checkbox"/>
72	<b>Absolute Peak Efficiency, Typical (%):</b>
>70	<b>Peak Efficiency Avg, Typical (%):</b>

## Regulatory Compliance

<a href="#">Compliant</a>	<b>RoHS 2015:</b>
<a href="#">View</a>	<b>Certificate of Conformance:</b>
<a href="#">Compliant</a>	<b>Reach 247:</b>
United States	<b>Country of Origin:</b>
Edmund Optics India Private Limited 267, Greystone Building, Second Floor, 6th Cross Rd, Binnamangala, Stage 1, Indiranagar, Bengaluru, Karnataka, India 560038 Phone: +91- 80-6845 0000	<b>Imported By:</b>

## Product Details

- Up to 80% Efficiency
- Bare Aluminum Coating
- Blaze Wavelength Options from 250 to 1600nm

Reflective Ruled Diffraction Gratings are optical components that separate (diffract) white (polychromatic) light into its component wavelengths. Theoretical resolving power is directly proportional to the number of grooves per millimeter ruled into the master grating. The blaze angle and groove spacing determine the specific wavelength and spectral region where the grating has its maximum efficiency. Reflective Ruled Diffraction Gratings are high quality, instrument-grade gratings that are ideal when high efficiency is the primary concern in devices such as monochromators and spectrometers. Each grating is fabricated from a highly accurate master grating that is copied many times. The duplication process is described below for replicated gratings.


Diffraction ruled gratings are used in a variety of monochromators for research, student and industry use. Almost all commercially available spectrophotometers (ultraviolet, visible, infrared, fluorescence, Raman, atomic absorption) utilize diffraction gratings to select specific wavelengths or scan over a wavelength interval. Generally, replicated ruled gratings should be used when high peak efficiency and throughput is required. Replicated holographic gratings should be used when minimum stray light is critical and high resolution is needed. Note: Damage thresholds for both ruled and holographic gratings are 350 milli-joules/cm<sup>2</sup> for pulsed lasers and 40 Watts/cm<sup>2</sup> for CW lasers.

**Handling Gratings:** Gratings require special handling, making them prone to fingerprints and aerosols. Gratings should only be handled by the edges. Before attempting to clean a grating, please [contact us](#).

### Manufacturing of Ruled Gratings

Ruled diffraction gratings are produced by ruling a series of closely spaced, straight parallel grooves into an optically flat aluminum coated substrate, known as the "master grating." Precise, interferometrically-controlled ruling engines utilize a very fine diamond tool to form a sawtooth-shaped groove profile at a given angle (commonly referred to as the blaze angle) on the surface of a prepared substrate. The replication process begins with the surface contour of a ruled master grating being vacuum deposition-coated with an extremely thin separation layer. An aluminum coating is then deposited on top of this separation layer. Then, an epoxy-coated flat glass substrate is placed on top of the layer-covered master, duplicating the grooved surface. The combination is cured and the process is finished when the replicated grating is separated from the master grating.

## Technical Information

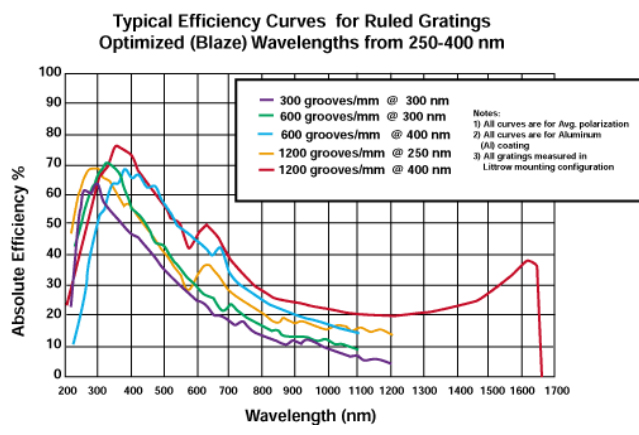

## Special Handling

These optics require special handling to avoid damage and ensure long-term performance. Proper handling, cleaning, and storage are essential to maintain optical quality. Explore our [Optics Cleaning Resources](#) for step-by-step guides and best practices. For personalized assistance, [Email us](#) or [Chat](#) with our technical support team.



Component Handling Tools