

[See all 20 Products in Family](#)

TECHSPEC® Fluorescence Filter Set for Cyan Fluorescence Protein (CFP)



Fluorescence Filter Sets



Stock #67-002 [CONTACT US](#)

1 ₹65,005

[ADD TO CART](#)

Volume Pricing	
Qty 1-9	₹65,005 each
Qty 10-25	₹57,220 each
Qty 26-49	₹50,992 each
Need More?	Request Quote

Product Downloads

General

Fluorescence Filter Kit **Type:**

CFP **Compatible Fluorophore:**

[#67-026](#) **Excitation Filter:**

Emission Filter:
#67-028

Dichroic Filter:
#67-078

Optical Properties

Substrate:
Fused Silica (Corning 7980)

Coating:
Hard Coated

Excitation Wavelength (nm):
426 - 450

Emission Wavelength (nm):
467 - 498

Dichroic Cut-On Wavelength (nm):
458.00

Wavelength Range (nm):
426 - 498

Regulatory Compliance

Certificate of Conformance:
[View](#)

Country of Origin:
United States

Imported By:
Edmund Optics India Private Limited

Product Details

The TECHSPEC® Fluorescence Filter Set for Cyan Fluorescence Protein (CFP) contains a dichroic filter, an emission filter, and an excitation filter designed for the Cyan Fluorescence Protein (CFP) spectrum. The emission filter has a wavelength range of 502 - 538nm to pass wavelengths within the CFP emission spectrum while the excitation filter has a wavelength range of 426 - 450nm to pass wavelengths within the CFP excitation/absorption spectrum.


- Excitation, Emission, and Dichroic Filters for Fluorescence Imaging
- >93% Transmission and OD 6 Blocking for Maximum Brightness and Contrast
- Hard Sputtered Coatings on Single Substrate
- [TECHSPEC Pre-Mounted Filter Cubes](#) Also Available

Our TECHSPEC® Fluorescence Filter Sets are designed to provide brilliant images when integrated into a fluorescence microscope. The high transmission of the bandpass filters guarantees the brightest possible images, while the deep blocking ensures no unwanted light gets to the detector. This ensures images will be brighter and will have much greater signal-to-noise than what is commonly available with soft coated, evaporated filters. TECHSPEC Fluorescence Filter Sets are ideal for fluorescence microscopy or fluorescence imaging applications and are available for common fluorescent dyes and fluorescent proteins. Additionally these filters are also used in polymerase chain reaction (PCR) diagnostic instruments.

Fluorescence filter sets include an Excitation and Emission filter from our [TECHSPEC® Fluorescence Bandpass Filters](#) and a Dichroic filter from our [TECHSPEC® Dichroic Filters](#). Excitation filters are placed within the illumination path of fluorescence microscopes to pass only wavelengths within the excitation range of the fluorophore under inspection. Emission filters are placed within the imaging path of a fluorescence microscope to pass only wavelengths within the emission range of the fluorophore. Dichroic filters are placed in between the excitation filter and the emission filter to reflect the excitation signal towards the fluorophore while transmitting the emission signal to the detector.

Filter sets are listed with the most common fluorophores. For a full list of compatible fluorophores with each set, view [Fluorophores and Optical Filters for Fluorescence Microscopy](#). TECHSPEC Fluorescence Filter Sets are available pre-mounted in [Olympus or Nikon Filter Cubes](#).

Technical Information

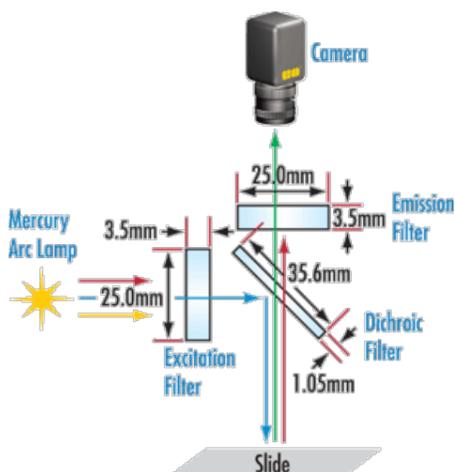

	Excitation and Emission Filters	Dichroic Filter
Size	25.0mm	25.2 x 35.6mm
Thickness	3.5mm	1.05mm
Tolerances		
Size	+0.0/-0.1	
Thickness	±0.1mm	
Clear Aperture	22.0mm	22.0 x 30.0mm
Transmitted Wavefront (RMS)	<λ/4	
Angle of Incidence	0.0° ±5.0°	45.0° ±5.0°

;