

[See all 20 Products in Family](#)

TECHSPEC® Fluorescence Filter Set for Qdot 525



Fluorescence Filter Sets



Stock **#86-374** [CONTACT US](#)

⊖ 1 ⊕ **MRP ₹95,846**

ⓘ Price inclusive of all taxes

ADD TO CART

| Volume Pricing | |
|----------------|-------------------------------|
| Qty 1-9 | ₹95,846 each |
| Qty 10-25 | ₹86,262 each |
| Qty 26-49 | ₹77,685 each |
| Need More? | Request Quote |

Product Downloads

General

Fluorescence Filter Kit **Type:**
Qdot® 525 **Compatible Fluorophore:**

#86-349

Excitation Filter:

#86-354

Emission Filter:

#86-332

Dichroic Filter:

Optical Properties

Fused Silica (Corning 7980)

Substrate: □

Hard Coated

Coating:

415 - 455

Excitation Wavelength (nm):

518 - 533

Emission Wavelength (nm):

510.00

Dichroic Cut-On Wavelength (nm):

415 - 533

Wavelength Range (nm):

Regulatory Compliance

[View](#)

Certificate of Conformance:

United States

Country of Origin:

Imported By:

Edmund Optics India Private Limited
267, Greystone Building, Second Floor,
6th Cross Rd, Binnamangala,
Stage 1, Indiranagar, Bengaluru,
Karnataka, India 560038
Phone: +91- 80-6845 0000

Product Details


- Excitation, Emission, and Dichroic Filters for Fluorescence Imaging
- >93% Transmission and OD 6 Blocking for Maximum Brightness and Contrast
- Hard Sputtered Coatings on Single Substrate
- [TECHSPEC Pre-Mounted Filter Cubes](#) Also Available

Our TECHSPEC® Fluorescence Filter Sets are designed to provide brilliant images when integrated into a fluorescence microscope. The high transmission of the bandpass filters guarantees the brightest possible images, while the deep blocking ensures no unwanted light gets to the detector. This ensures images will be brighter and will have much greater signal-to-noise than what is commonly available with soft coated, evaporated filters. TECHSPEC Fluorescence Filter Sets are ideal for fluorescence microscopy or fluorescence imaging applications and are available for common fluorescent dyes and fluorescent proteins. Additionally these filters are also used in polymerase chain reaction (PCR) diagnostic instruments.

Fluorescence filter sets include an Excitation and Emission filter from our [TECHSPEC® Fluorescence Bandpass Filters](#) and a Dichroic filter from our [TECHSPEC® Dichroic Filters](#). Excitation filters are placed within the illumination path of fluorescence microscopes to pass only wavelengths within the excitation range of the fluorophore under inspection. Emission filters are placed within the imaging path of a fluorescence microscope to pass only wavelengths within the emission range of the fluorophore. Dichroic filters are placed in between the excitation filter and the emission filter to reflect the excitation signal towards the fluorophore while transmitting the emission signal to the detector.

Filter sets are listed with the most common fluorophores. For a full list of compatible fluorophores with each set, view [Fluorophores and Optical Filters for Fluorescence Microscopy](#). TECHSPEC Fluorescence Filter Sets are available pre-mounted in [Olympus or Nikon Filter Cubes](#).

Technical Information

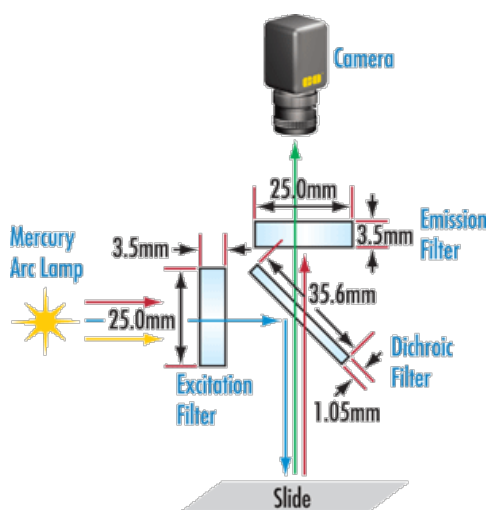

| | Excitation and Emission Filters | Dichroic Filter |
|------------|---------------------------------|-----------------|
| Size | 25.0mm | 25.2 x 35.6mm |
| Thickness | 3.5mm | 1.05mm |
| Tolerances | | |
| Size | +0.0/-0.1 | |
| Thickness | ±0.1mm | |

| | | |
|-----------------------------|---|------------------|
| Clear Aperture | 22.0mm | 22.0 x 30.0mm |
| Transmitted Wavefront (RMS) | $\lambda/4$ | |
| Angle of Incidence | 0.0° \pm 5.0° | 45.0° \pm 5.0° |
| Compatible Microscopes | Most common Nikon, Olympus, and Zeiss Microscopes | |
