


Hydrogen Cyanide Fiber-Coupled Gas Cell, 16.5 cm Path Length, SC/PC

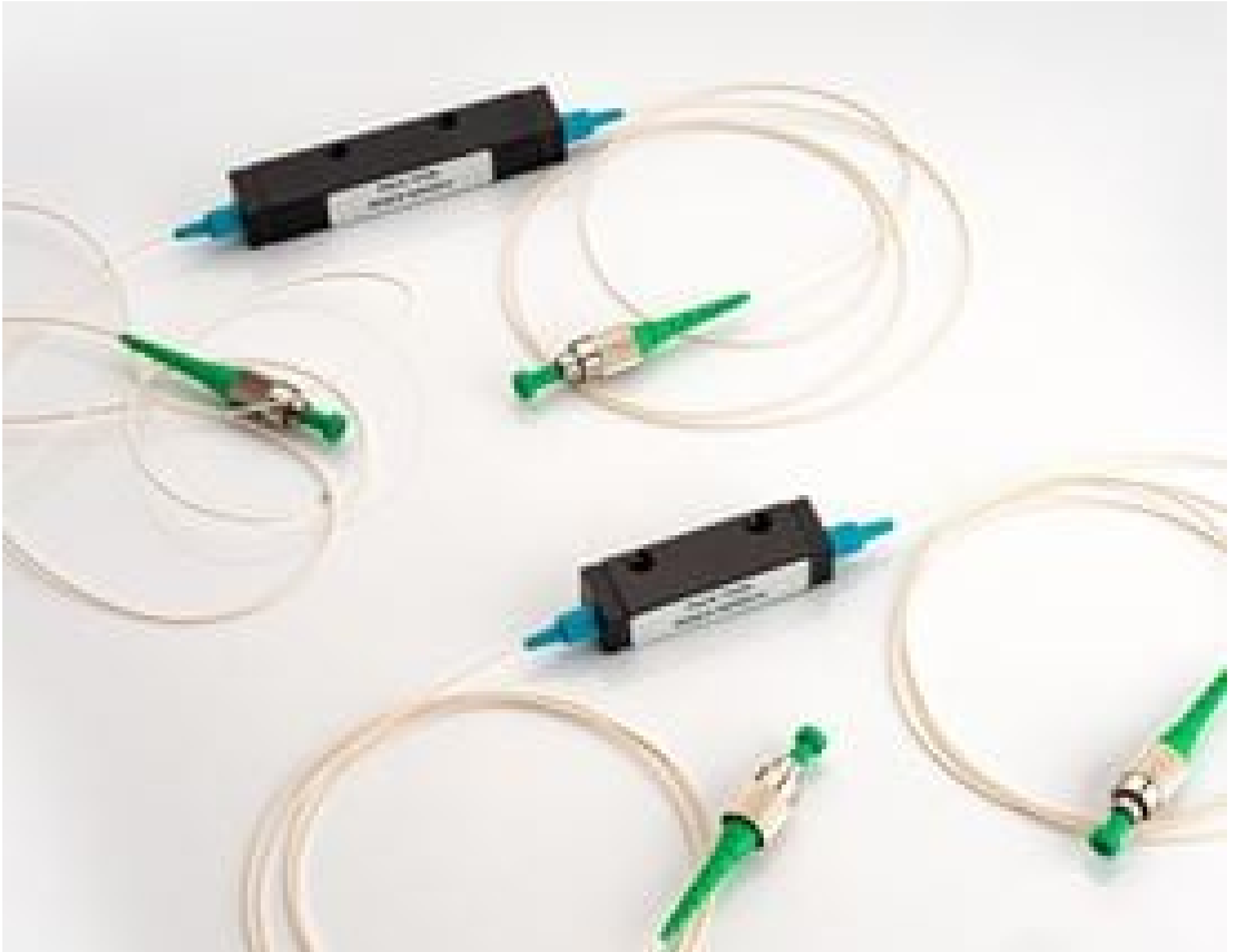

Stock #72-203 **1 In Stock**

- 1 + MRP ₹1,08,855

● Price inclusive of all taxes

ADD TO CART

Volume Pricing	
Qty 1-4	₹1,08,855 each
Qty 5-9	₹97,970 each
Need More?	Request Quote

Product Downloads

16.5 Path Length (cm):

SC/PC Fiber Connector Type:

General

Hydrogen Cyanide (H¹³ C¹⁴ N) Type:

HCN-13-H(16.5)-25-SCPC Model Number:

Optical Properties

Wavelength Range (nm):

1525 - 1565

Transmission (%):

>50

Environmental & Durability Factors

Operating Temperature (°C):

0 to +70

Regulatory Compliance

Certificate of Conformance:

[View](#)

Country of Origin:

United States

Imported By:

Edmund Optics India Private Limited
267, Greystone Building, Second Floor,
6th Cross Rd, Binnamangala,
Stage 1, Indiranagar, Bengaluru,
Karnataka, India 560038
Phone: +91- 80-6845 0000

Product Details

- Hydrogen Cyanide and Acetylene Gasses Available for Wavelength Coverage of 1510 to 1565nm
- FC and SC Fiber Connection Options with a Variety of Path Lengths Available
- NIST Traceable Reference Cells

Wavelength References Fiber-Coupled Gas Cells are FC/APC, SC/APC, FC/PC, or SC/PC connected fiber coupled, gas filled precision filters whose absorption wavelengths depend on specific molecular energy level transitions that may be used as wavelength standards. Hermetically sealed to maintain >10 year lifetime, these gas cells feature metal housings, wedged windows, and coated optics for minimum interference artifacts and can be easily integrated into existing benchtop systems. Wavelength References Fiber-Coupled Reference Cells are available with a variety of path lengths and pressures which meet the requirements of NIST Standard Reference Material® (SRMs) 2517a, 2519, or 2519a. Short path lengths are recommended for measuring gasses at high concentration while longer path lengths enable more sensitive measurements. These reference cells are ideal for spectroscopy, wavelength/frequency locking, laser calibration, and optical gas sensing systems. Hydrogen Cyanide ($H^{13}C^{14}N$) has been identified by national standards bodies as the primary wavelength reference in the C-band (1530 – 1565nm) while Acetylene ($^{12}C_2H_2$) is recognized as a primary wavelength reference in the 1510 to 1540nm band.