

☐ - 1 ☐ + ₹21,750

Volume Pricing	
Qty 1+	₹21,750 each
Need More?	Request Quote

General	
1	General
2	General
3	General
4	General
5	General
6	General
7	General
8	General
9	General
10	General
11	General
12	General
13	General
14	General
15	General
16	General
17	General
18	General
19	General
20	General
21	General
22	General
23	General
24	General
25	General
26	General
27	General
28	General
29	General
30	General
31	General
32	General
33	General
34	General
35	General
36	General
37	General
38	General
39	General
40	General
41	General
42	General
43	General
44	General
45	General
46	General
47	General
48	General
49	General
50	General
51	General
52	General
53	General
54	General
55	General
56	General
57	General
58	General
59	General
60	General
61	General
62	General
63	General
64	General
65	General
66	General
67	General
68	General
69	General
70	General
71	General
72	General
73	General
74	General
75	General
76	General
77	General
78	General
79	General
80	General
81	General
82	General
83	General
84	General
85	General
86	General
87	General
88	General
89	General
90	General
91	General
92	General
93	General
94	General
95	General
96	General
97	General
98	General
99	General
100	General

Type: Camera Accessory

Certificate of Conformance:

Edmund Optics India Private Limited

Product Details

- Long Exposure CMOS (Over 60 Minutes)
- Feature Extended Dynamic Range (EDR)
- Dual Imaging Modes for Flexibility
- Scientific Cameras for High Sensitivity Microscopy Applications




Teledyne Photometrics Retiga E USB Cameras feature ultra-low dark current (<0.001 electrons/pixel/second), enabling exposure times of over an hour for low-light and long-duration imaging. The cameras feature EDR mode, which combines a large full-well capacity with low read noise to capture both intense and weak signals in the same image. The Retiga E series offers dual imaging modes; a long exposure mode enables ultra-sensitive, low-light imaging, and a fast capture mode for high-speed acquisition. The dual modes provide flexible performance for both static and dynamic scientific applications. Teledyne Photometrics Retiga E USB Cameras are ideal for time-lapse imaging in cell colony growth, fluorescence microscopy, and bioluminescence/chemiluminescence imaging. These cameras include Teledyne's proprietary-designed software platforms, Beacon and PVCAM, for optimizing camera performance and ease of system integration.

Long Exposure Imaging


Thanks to major breakthroughs in thermal noise control, Retiga E CMOS cameras have an unprecedented ultra-low dark current of <0.001 electrons per pixel per second, allowing for exposure times of over an hour.

Stretching exposure times from seconds to minutes to hours in this manner results in comfortable detection of even the weakest signals, ideal for imaging luminescence.

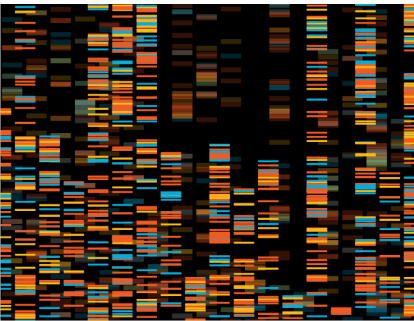

Extended Dynamic Range

As well as Speed and Long Exposure modes, Retiga E cameras feature Extended Dynamic Range (EDR) mode, which combines a large full well capacity with a low read noise to result in ~2.5x more dynamic range than other modes, easily capturing both intense and weak signals in the same image.

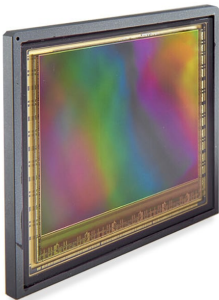


Small Pixel Imaging

The Retiga E20 features our smallest-ever pixel at 2.4 μm , which results in high-resolution imaging at even the lowest magnifications. Combined with a 20-Megapixel array (5472 x 3648) the Retiga E20 can pick out even the smallest details from large samples in a high-throughput manner.

With low noise contribution and a high peak QE, the Retiga E20 will deliver high quality images of your samples even in low-light conditions.

High Speed Mode

As well as the 'Long Exposure' mode, Retiga E cameras feature a 'Fast Capture' mode that operates at 100 fps for the Retiga E7, across the full 7 megapixel sensor. This flexibility allows the Retiga E7 to capture fast, dynamic signals, or to switch to long exposure in order to service a wide range of imaging applications.

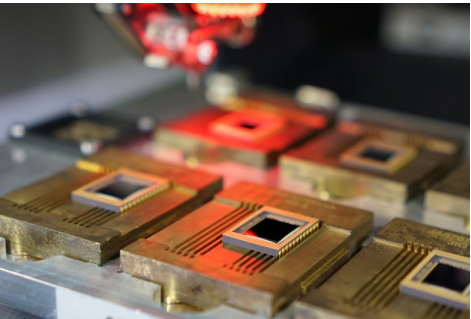

OEM Integration

Interfacing is simple with the Retiga E Family by design. With their compact format and single cable interface, hardware integration is straightforward. The addition of onboard triggering, Windows/Linux support, and our dedicated OEM team make integrating and customising the Retiga E7 or Retiga E20 a painless process.