

Coherent® Lasercam™ USB Laser Beam Profiler 1282870 | 2/3" Format

See More by [Coherent®](#)



Stock #33-157 **1 In Stock**

MRP ₹6,14,898

Price inclusive of all taxes

ADD TO CART

Volume Pricing	
Qty 1+	₹6,14,898 each
Need More?	Request Quote

Product Downloads

Physical & Mechanical Properties

Dimensions (mm):
68.1 x 79.3 x 40.9 (with LDFP)

Weight (g):
110 (without cable)

Dimensional Accuracy:
±1% (typical), ±5% max (over entire spectral and dimensional range)

Optical Properties

Beam Diameter, 1/e² (mm):

Recommended: 0.20mm min, 6.0mm max

Spectral Range:
300 - 1100nm (400 - 1100nm with LDFP)

CW Saturation:
5 mW/cm² (at 632.8nm with LDFP)

Pulsed Saturation:
0.15 mJ/cm² (at 1.06µm with LDFP)

Peak Intensity:
Recommended: 75-95% of camera saturation

Damage Threshold, By Design:
32 mJ/cm² @ 1.06µm without Low Distortion
Faceplate

Sensor

Sensing Area, H x V (mm):
8.3 x 6.6

Gamma:
1.00

Sensor Format:
2/3"

Frame Rate:
27Hz (Live Mode), 10Hz (with calculations)

Exposure Time:
Fixed at 10ms

Electrical

Signal to Noise S/N Ratio (dB):
>60

Gain (dB):
Gain is factory set for optimum linear dynamic range, not user adjustable

Peak Noise (nW/cm²):
24 (at 632.8nm)

Pulse Trigger:
TTL, rising or falling edge

Hardware & Interface Connectivity

Operating System:
Windows®

Threading & Mounting

Mount:
C-Mount

Environmental & Durability Factors

Operating Temperature (°C):
-20 to 60

Regulatory Compliance

RoHS 2015:
[Exempt](#)

Reach 224:
[Contains SVHC\(s\)](#)

Certificate of Conformance:
[View](#)

Country of Origin:
United States


Imported By:
Edmund Optics India Private Limited
267, Greystone Building, Second Floor,
6th Cross Rd, Binnamangala,
Stage 1, Indiranagar, Bengaluru,
Karnataka, India 560038
Phone: +91- 80-6845 0000

Product Details

- 12 and 14-bit Digital USB 2.0 Interface Options
- High Sensitivity and Dynamic Range
- Intuitive BeamView™ Software Included

The Coherent® Lasercam™ Beam Profiler features excellent signal-to-noise ratio and linear response for accurate pulsed as well as CW laser beam dimension and uniformity measurements. The new BeamView™ 4.4 interface software features TCP/IP control and NI LabVIEW™ library suite, enabling efficient and smooth integration of beam profiling into any application. Examples of analysis functions that can be performed on the acquired beam images include: beam centroid location, beam peak intensity position, pointing stability, total relative power/energy in beam, peak power/energy density of beam, beam divergence, ellipticity, beam intensity uniformity, Gaussian fit, beam diameter/width based on second moments or user selectable percentage of peak/total energy.

Technical Information

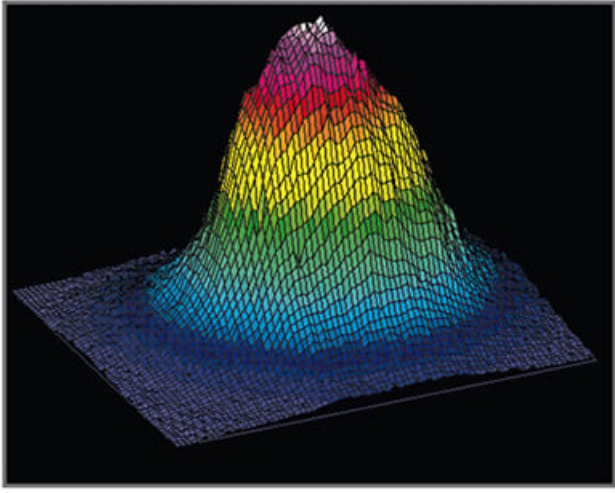

Intuitive Software Interface

