

Coherent® Lasercam™ USB Laser Beam Profiler 1282870 | 2/3" Format

See More by [Coherent®](#)



Stock **#33-157** **2 In Stock**

-

1

+

₹5,38,200

ADD TO CART

Volume Pricing	
Qty 1+	₹5,38,200 each
Need More?	Request Quote

Product Downloads

Physical & Mechanical Properties

Dimensions (mm):	
68.1 x 79.3 x 40.9 (with LDFFP)	
Weight (g):	
110 (without cable)	
Dimensional Accuracy:	
±1% (typical), ±5% max (over entire spectral and dimensional range)	
Optical Properties	
Beam Diameter, 1/e ² (mm):	
Recommended: 0.20mm min, 6.0mm max	

300 - 1100nm (400 - 1100nm with LDFFP)	Spectral Range:
5 mW/cm ² (at 632.8nm with LDFFP)	CW Saturation:
0.15 mJ/cm ² (at 1.06µm with LDFFP)	Pulsed Saturation:
Recommended: 75-95% of camera saturation	Peak Intensity:
32 mJ/cm ² @ 1.06µm without Low Distortion Faceplate	Damage Threshold, By Design: ☐

Sensor	
8.3 x 6.6	Sensing Area, H x V (mm):
1.00	Gamma:
2/3"	Sensor Format:
27Hz (Live Mode), 10Hz (with calculations)	Frame Rate:
Fixed at 10ms	Exposure Time:

Electrical	
>60	Signal to Noise S/N Ratio (dB):
Gain is factory set for optimum linear dynamic range, not user adjustable	Gain (dB):
24 (at 632.8nm)	Peak Noise (nW/cm²):
TTL, rising or falling edge	Pulse Trigger:

Hardware & Interface Connectivity	
Windows®	Operating System:

Threading & Mounting	
C-Mount	Mount:

Environmental & Durability Factors	
-20 to 60	Operating Temperature (°C):


Regulatory Compliance	
Exempt	RoHS 2015:
Contains SVHC(s)	Reach 224:
View	Certificate of Conformance:
United States	Country of Origin:
Edmund Optics India Private Limited	Imported By:

Product Details

- 12 and 14-bit Digital USB 2.0 Interface Options
- High Sensitivity and Dynamic Range
- Intuitive BeamView™ Software Included

The Coherent® Lasercam™ Beam Profiler features excellent signal-to-noise ratio and linear response for accurate pulsed as well as CW laser beam dimension and uniformity measurements. The new BeamView™ 4.4 interface software features TCP/IP control and NI LabVIEW™ library suite, enabling efficient and smooth integration of beam profiling into any application. Examples of analysis functions that can be performed on the acquired beam images include: beam centroid location, beam peak intensity position, pointing stability, total relative power/energy in beam, peak power/energy density of beam, beam divergence, ellipticity, beam intensity uniformity, Gaussian fit, beam diameter/width based on second moments or user selectable percentage of peak/total energy.

Technical Information

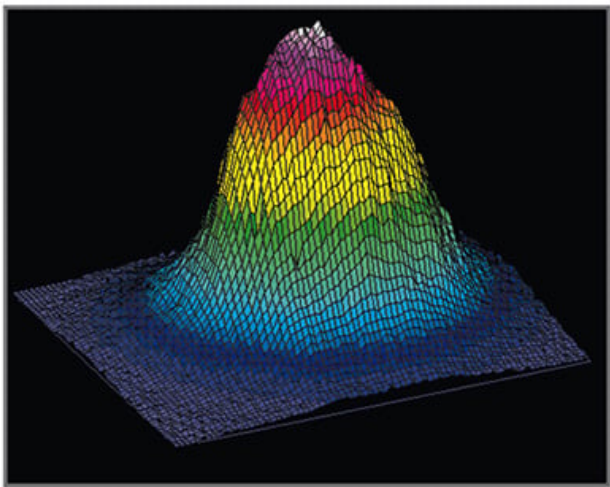

Intuitive Software Interface

